



The scentillating choice for N<sub>2</sub>O delivery... because patient comfort is your top priority.



## Personal Inhaler Plus+™ Nasal Hoods

- Reduces risk of cross-contamination
- Relaxes patients, optimizing operator control
- Saves valuable staff time
- 8 soothing scents/3 comfortable sizes
- Made in USA for maximum quality standards
- Non-latex material minimizes potential for allergic reactions
- Thermoplastic elastomer reduces concern for nasal hood collapse

## Accutron Scavenging Circuit

- ▶ Easy to-read vacuum-flow control gauge
- ▶ Low-profile scavenging assembly expands clinical field of operation
- ▶ Clear exhalation valve hub permits respiratory monitoring
- ▶ All components latex-free
- ▶ Convenient X-ray compatible circuitry
- ▶ Easy-to-clean, autoclavable circuitry
- ▶ Flexible system accommodates all brands of flowmeters
- ▶ Bright white finish provides an updated clinical appearance
- ▶ Made in U.S.A. for maximum quality standards



# PIP+™ Nasal Hoods

Accutron's Personal Inhaler Plus+™ single-patient-use nasal hood gives your patients reliable, comfortable, high-quality N<sub>2</sub>O delivery.

As the only colored and scented nasal hood manufactured in the USA, Accutron's PIP+ sets the quality standard. The PIP+ is made of thermoplastic elastomer to minimize concern for nasal hood collapse and is certified latex-free, reducing the possibility of allergic reactions.



Packages of 24 brightly colored nasal hoods are offered in eight different scents (plus unscented gray). Available in three sizes—large, medium, small—the PIP+ provides patients with a comfortable fit and a desirable clinical seal to prevent N<sub>2</sub>O leakage.

Take special care of your patients, your staff and yourself when you offer your patients the Accutron Personal Inhaler Plus+.



The original flavor-scented nasal hood that pleases patients and eases your procedures.

## Order PIP+™ Nasal Hoods in these scents and sizes

### Small - 24 Per Package

- 33017-17 Groovy Grape
- 33017-12 Birthday Bubblegum
- 33017-11 Sassy Strawberry
- 33017-10 Outlaw Orange
- 33017-9 Unscented Whisper Gray
- 33017-1 Variety Pack 1 (Orange, Strawberry & Bubblegum)

Nasal hoods are individually wrapped and packaged in 24-count containers. To maintain scent integrity, Variety Pack same-scent hoods are packaged in individual zip-locked bags.

### Medium - 24 Per Package

- 33016-16 Fresh Mint
- 33016-15 Perfect Peach
- 33016-14 French Vanilla
- 33016-13 Piña Colada
- 33016-12 Birthday Bubblegum
- 33016-11 Sassy Strawberry
- 33016-10 Outlaw Orange
- 33016-9 Unscented Whisper Gray
- 33016-1 Variety Pack 1 (Orange, Strawberry & Bubblegum)
- 33016-2 Variety Pack 2 (Piña Colada, Vanilla & Orange)
- 33016-3 Variety Pack 3 (Mint, Peach & Vanilla)

### Large - 24 Per Package


- 33015-16 Fresh Mint
- 33015-15 Perfect Peach
- 33015-14 French Vanilla
- 33015-9 Unscented Whisper Gray
- 33015-3 Variety Pack 3 (Mint, Peach & Vanilla)

**24 nasal hoods per package**

See back page for Component Products & Adaptor Easy Order Guide.



## Scavenging Circuit I (33010)

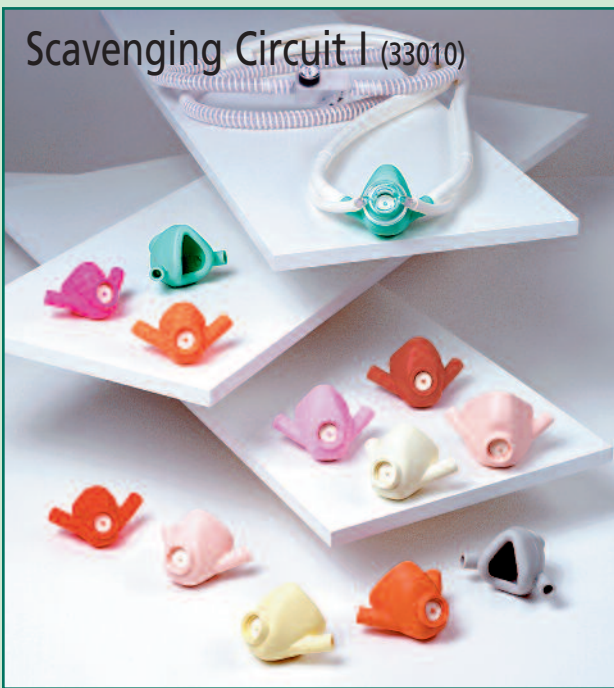

# Scavenging Circuit

Accutron's **Scavenging Circuit** reduces residual nitrous oxide in treatment rooms and surrounding areas, providing patients and staff with the safest environment possible.

X-ray compatible and autoclavable, the system consists of a clear exhalation valve hub, a vacuum control valve assembly (33010 & 32007 only) and a continuum of tubing. The transparent hub is unique in that it permits visual monitoring of patient respiration during dental procedures. In addition, the hub's low-profile design offers an expanded field of operation, a definite plus since treatment space is often in short supply. The circuit's vacuum control assembly meets NIOSH and ADA proposed scavenging flow rate recommendations. And its gauge is easy


to read and highly visible.

For your  
protection.

With its bright white tubing, the circuit creates a fresh clinical appearance as it safely and effectively delivers and removes nitrous oxide. Add an Accutron PIP+ single-patient-use nasal hood to ensure proper facial seal and complete the system.

**Accutron Scavenging Circuit**—the surest way to minimize patient anxiety and secure a safe and productive work environment for doctor and staff. Latex-free components.

## Scavenging Circuit II (32007)

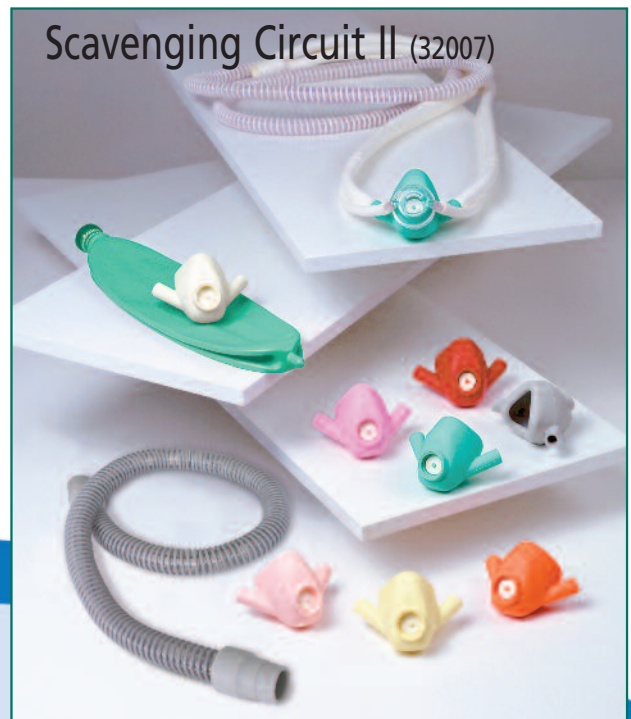

### Order No. Description

- |       |   |
|-------|---|
| 33010 | <b>Scavenging Circuit I:</b> Scavenging Circuit with Vacuum Control Valve and complimentary sampling of single-use PIP+ Nasal Hoods (nasal hoods included are representative of all sizes available)  |
| 32007 | <b>Scavenging Circuit II:</b> Scavenging Circuit with Vacuum Control Valve, Corrugated Tube, Reservoir Bag and complimentary sampling of single-use PIP+ Nasal Hoods (nasal hoods included are representative of all sizes available)           |
| 33008 | <b>Scavenging Circuit I for Remote Flow System and Accu-Vac:</b> Scavenging Circuit and complimentary sampling of single-use PIP+ Nasal Hoods (nasal hoods included are representative of all sizes available)                                  |
| 33009 | <b>Scavenging Circuit II for Remote Flow System and Accu-Vac:</b> Scavenging Circuit, Corrugated Tube, Reservoir Bag and complimentary sampling of single-use PIP+ Nasal Hoods (nasal hoods included are representative of all sizes available) |

See back page for Component Products.

# Component Products

## Adaptors/Reducers – PIP+™ & Multi-Use Nasal Hoods to Scavengers:

- 35000 Porter Brown Adaptor
- 35100 Porter Clean Air Adaptor
- 35325 MatrX Spiral Tubing Adaptor
- 35350 MatrX Smooth Tubing Adaptor
- 35400 McKesson Adaptor
- 35600 Non-Scavenging System Tube Size Reducing Kit
- 20305 Side Tube Reducers (2/Pkg)
- 26250 Cone Ring (Matrx)

## Scavenging Circuits with Multi-Use Nasal Hoods:

- 33011 Scavenger with Large Multi-Use Hood
- 33012 Scavenger with Medium Multi-Use Hood
- 33013 Scavenger with Small Multi-Use Hood



## Multi-Use Nasal Hoods:


- 33018 Multi-Use Large Hood
- 33019 Multi-Use Medium Hood
- 33020 Multi-Use Small Hood

## Scavenging Circuit Components:

- 21178 Double Connectors for Scavenging System (2/Pkg)
- 21179 Single Connectors for Nasal Hood (2/Pkg)
- 21245 Large Y Connector for Scavenging System
- 21246 Small Y Connector for Scavenging System
- 22099 Scavenging Hub (clear disk)
- 22122 Vacuum Gauge
- 25591 Spiral Vacuum Tube – 3.5'
- 25592 Spiral Vacuum Tube – 1.5'
- 25593 5" Large Tube
- 25594 7 1/4" Small Tube
- 26253 Tube Slide Double
- 26254 Single Vacuum Connector
- 26255 90° Elbow (creates 90° angle off flowmeter for attachment of scavenging circuit)



## Scavenging Circuit Parts Diagram

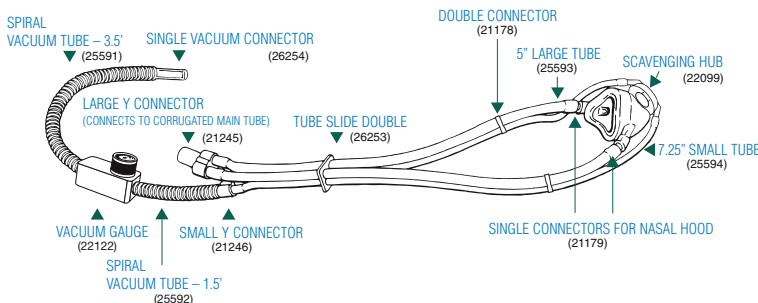

**Note:** If using RFS & Accu-Vac Scavenging Circuit, parts 25591, 22122 and 25592 are replaced with a 6' Spiral Vacuum Tube (part 26196).

## Rubber Goods:

- 22137 6-foot Corrugated Tube
- 33014 3-foot Corrugated Tube
- 33025 Reservoir Bag – 3 liter
- 33026 Reservoir Bag – 2 liter

# Adaptor Easy Order Guide

### Porter Brown Scavenging System - 35000

Has 4 tubes - 2 ports in each side, hood within hood bearing the name "Brown."

### Porter Clean Air Scavenging System - 35100

Has tube inside tube and wings or vacuum sweeps to draw off excess gas.

### Matrx Scavenging System - 35350

Has one smooth tube on each side.

### Matrx Scavenging System - 35325

Has spiral, white tubing.

### McKesson Scavenging System - 35400

Has clear hood and 1 tube on each side.

### Non-Scavenging System Tube Size Reducing Kit - 35600

Reduces tubing to proper size for easy connection.

### Accutron Scavenging System does not require an adaptor.

